



## Policy & Procedure (P & P)

Policy Title :

### Screening of Donated Blood

Department	Index No.	Scope
Laboratory & Blood Bank	LAB-074	All Blood Bank and Lab Staff
Issue Date	Revision NO	Effective Date
1432/06/04	4	1440/07/23
Review Due Date	Related Standard NO.	Page Number#
1442/07/23	CBAHI ( LB.51 )	4

#### 01. Policy:

- 01.1. All blood units are screened to the Transfusion Transmitted Disease markers according to the AABB policy by ELISA or chemiluminescence assay AND Nucleic Acid Amplification Test: NAT tests.

#### 02. Definition :

N/A

#### 03. Purpose :

- 03.1. To provide safe blood free from Transfusion Transmissible Disease (TTD) according to the AABB policy.

#### 04. Procedure :

- The blood bank technician, every morning, writes the serial numbers of the Blood bank serology report with signature and stamp then sends the samples to the Serology department where serological tests are performed by chemiluminescence assay (Architect) and Nucleic Acid Technique "NAT".
- The serology department performs the following tests to each donated unit:
  1. HBs Ag.
  2. Anti-HBc
  3. Anti-HBs Ab for all Anti-HBc positive samples
  4. Anti-HCV





وزارة الصحة  
Ministry of Health  
مستشفى القنفذة العام

5. Anti-HTLV I, II
6. Serological test for syphilis (RPR).
7. HIV I/ II combo (anti HIV I/ II antibodies and P24 Antigene)
8. HIV RNA
9. HBV DNA
10. HCV RNA
11. Malaria thick film stained by Giemsa stain.

**Note:**

1. Donor's sample is collected in 6 ml plain tube to perform the chemiluminescence assay.
2. Donor's sample is collected in 7ml EDTA tube to perform NAT Test for (HIV/HBV/HCV) for all blood donors.
3. All blood component units are stored in a separate secured proper storage places: a separate refrigerator labeled UNSCREENED BLOOD until all TTD test results released then the bags will be shifted to another refrigerator for SCREENED BLOOD or discarded.
4. **Interpretation of Tests:**
  - 4.1. Results are received from serology department as reactive (+) or non-reactive (NR) for the tests.
  - 4.2. Two blood bank staff Compare the donor's results by both chemiluminescence assay and NAT tests.
  - 4.3. Blood is discarded if either chemiluminescence assay or NAT tests is positive.
  - 4.4. All blood component units with positive screen TTD test results are repeated in duplicate one sample from blood unit and another sample from the test tube before the release of the whole batch of donation, to solve any problems about labeling.
  - 4.5. All repeat TTD positive screen test result must be confirmed by NEUTRALISATION for HBV or Western Blot for (HIV I, II and HTLV I, II) and RIBA for HCV and PCR for NAT tests which was approved by FDA or CE with blood sample from the donor.
  - 4.6. All units positive for HBs Ag, HCV Ab, HIV I, II Ab, HIV P24, HTLV I, II, Syphilis, HIV RNA, HCV RNA, HBV RNA or malaria are discarded.
  - 4.7. Anti HBs test should be performed to all anti-HBc positive samples:



- if antiHBs result is negative discard blood units
  - if anti-HBs result is positive and the antibody titer is above 1000 mIU/ml keep the unit and use it due to presence of Natural Immunity in this donor for HBV regarding that DNA-NAT test is negative.
  - if anti-HBs result is positive and the antibody titer is between 100 mIU/ml and 1000 mIU/ml keep the units of platelets and discard the units of Fresh Frozen Plasma FFP and keep the unit of PRBCs in quarantine and use it only in emergencies.
- 4.8. Donors with positive TTD tests are referred to the Infection Control Department for follow up.
- 4.9. All blood donors with positive screening or confirmatory results for HIV I, II and HTLV I, II, HBV and HCV are DEFERRED PERMANENTLY.
- 4.10. All blood donors with positive MALARIA test are DEFERRED FOR THREE YEARS from the date of treatment and cessation of symptoms.
- 4.11. All blood donors with SYPHILIS test or history of Gonorrhea are DEFERRED FOR ONE YEAR from the date of treatment and cessation of symptoms.
- 4.12. All blood component units with negative screen TTD test results are rechecked by TWO QUALIFIED PERSONS from blood bank (BB supervisor or technician and physician) and should sign the result before labeling.
- 4.13. All blood component units with rechecked negative screen TTD test results are labeled with a special label sticker with the mention of NON-REACTIVE for all TTD then are stored in separate secured proper storage refrigerator labeled SCREENED BLOOD after release of all results to be ready for transfusion.
- 4.14. All blood component unit's storage places (Fridges, Freezers and platelet incubators) are labeled: SCREENED BLOOD for the screened Blood and UNSCREENED BLOOD for unscreened blood
- 4.15. On release of any Non-Reactive Blood Component units the technician will check again for TTD test results reports and check if it is non-reactive before release.





وزارة الصحة  
Ministry of Health  
مستشفى القنفذة العام

**05. Responsibilities :**

05.1. All Blood Bank Staff, Serology Staff and Hematology Staff of Al-Qunfudah General Hospital

**06. Equipment & Forms**

06.1. Blood Bank Report Form for serology results.

06.2. Blood Bank Report Form for malaria parasite.

**07. Attachment :**

N/A

**08. Reference**

08.1. The Technical Manual of the American Association of Blood Banks.

**Preparation , Reviewing & Approval Box**

	NAME	POSITION	SIGN & STAMP	DATE
Prepared By	Dr RAJA NACER SASSI	Head of Blood Bank		
Reviewed By	Mr. ABDULHADI ASHIRI	Lab & B. Bank HOD		
Document Reviewed By	Ms. SADIHA ALMAHMOUDI	TQM Director		12/10/20
Reviewed By	Dr. AGEEL ALGANIMI	Medical Director		
Approved By	Dr. ABDULLAH ALJABRI	Hospital Director		12/11/20








وزارة الصحة  
Ministry of Health

مستشفى القنفذة العام  
Al Qunfudah General Hospital

**Appendix A: (Donor screening flow chart ) :**

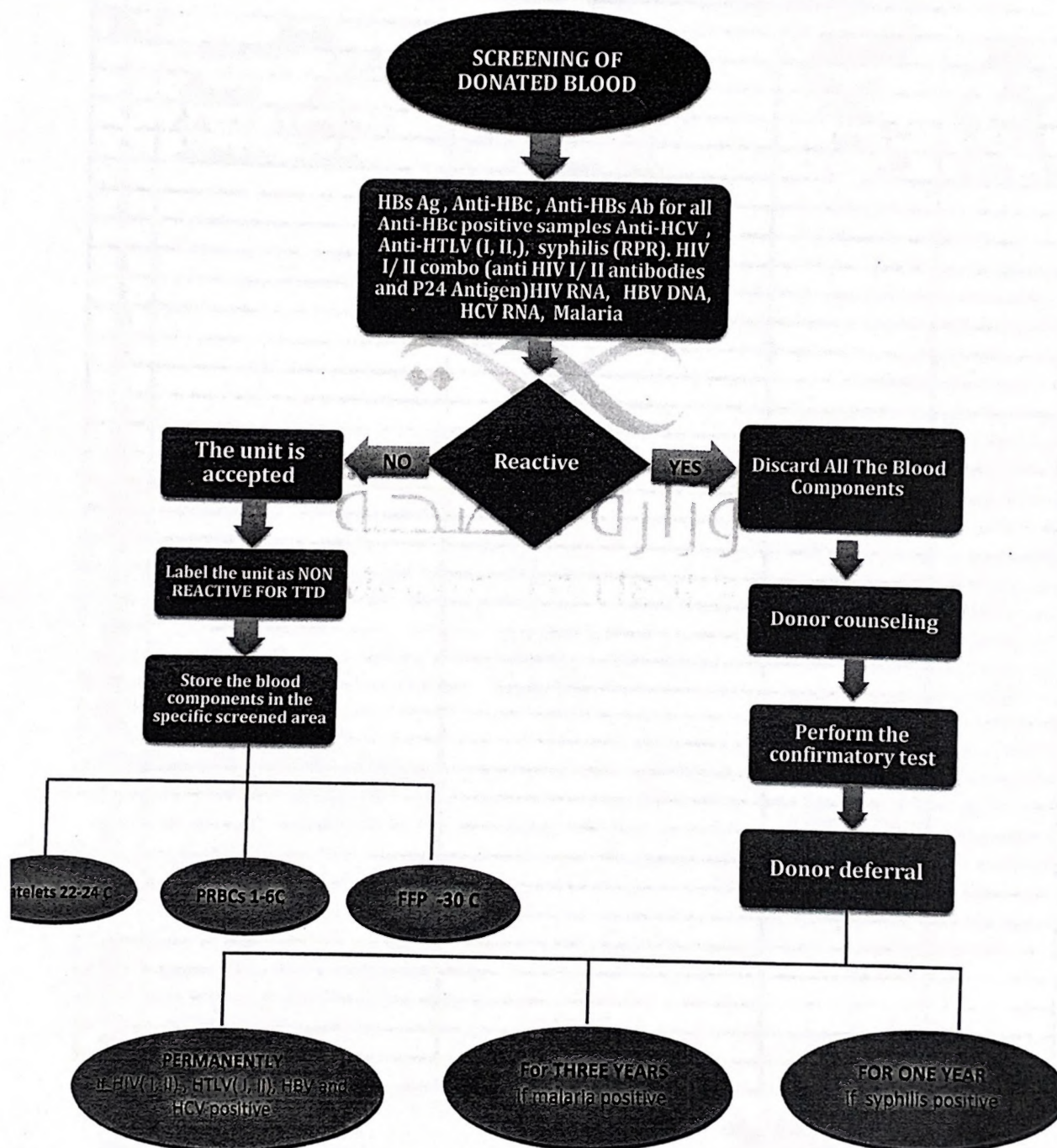





# AL QUNFUDAH GENERAL HOSPITAL

Department: Blood Bank

Signature Sheet For Policy And Procedure				
Policy And Procedure Title		Policy And Procedure No.		
Screening of Donated Blood		LAB-074		
No.	Name	Title	Signature	Date
1	Monemeh Alshamir	Lab Tech		12-6-19
2	Ahmad Ali	Lab Tech		12.6.19
3	Ahmad Alqabani	Lab Tech		12/6/19
4	Hassan Zahr	Lab Tech		12-6-19
5				
6				
7				
8				
9				
10				
11				
12				
13				
14				
15				
16				
17				
18				
19				
20				
21				
22				
23				
24				
25				
26				
27				
28				
29				
30				
31				
32				
33				
34				
35				
36				